



ACCENT<sup>®</sup>XL  
**PERIORBITAL**  
 H A N D P I E C E

The **NEW** PeriOrbital handpiece uniquely designed for UniPolar<sup>™</sup> RF-based treatments around the eyes

- ▶ **Improvement of fine lines and wrinkles in the periOrbital area**
- ▶ **Crow's feet** ▶ **Infraorbital area**

Based on UniPolar RF-based technology, the PeriOrbital handpiece provides the ideal tool for non-volumetric heating indicated for wrinkles around the eyes. With its perfect size and configuration for the treatment of skin in the infraorbital area and crow's feet, the unique mushroom-shaped tip makes Accent<sup>XL</sup> treatments on the face even easier. Try it and you'll see the difference!

▶ **IN-Motion Tip**

Ergonomically contoured IN-Motion tip allows treatment around the eyes in under 10 minutes.



**Innovative.** The PeriOrbital handpiece uses IN-Motion<sup>™</sup> technology to deliver UniPolar RF energy in the hard-to-treat areas around the eyes.

**Fast.** Ergonomically contoured IN-Motion tip allows treatment around the eyes in under 10 minutes.

**Affordable.** The PeriOrbital module is a simple upgrade to the Accent<sup>XL</sup>, no separate system is needed.

**Ease of Use.** Lightweight, ergonomic handpiece with a contoured tip makes treating around the eyes easy. No consumables required.

# Accent<sup>XL</sup>: The Most Advanced RF Multi-Application System from Alma Lasers




## Dual-layer Thermotherapy

The Accent<sup>XL</sup> is a unique, upgradeable, multi-application, Radio Frequency (RF) based platform for non-invasive skin tightening, body contouring, treatment of cellulite, and **now skin rejuvenation**. The BiPolar and UniPolar applications of the Accent<sup>XL</sup> work by gradually heating the dermal and subdermal tissue to smooth, tighten and re-contour both superficially **and** volumetrically.

## Technical Specification of PeriOrbital Module

<b>Energy Source</b>	RF Energy
<b>Frequency</b>	40.68 MHz
<b>Frequency Stability</b>	+/-5%
<b>Maximal RF-Power</b>	90 Watts
<b>Operative RF-Power</b>	40-60 Watts
<b>Pulse Repetition Rate</b>	15 kHz
<b>Material of Tip</b>	Aluminum coated with Alumina
<b>Contactable Area Diameter</b>	12 mm

The handpiece operates at two distinct cycle durations for a total of 20 seconds/cycle:

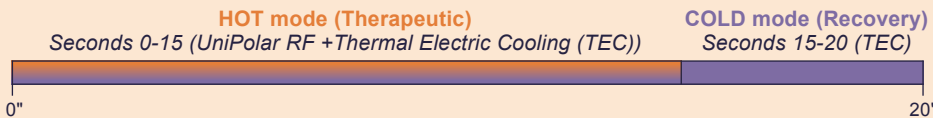

During the first 15 seconds, the UniPolar RF and TEC are active at the same time (HOT + COLD). Five seconds of TEC alone (COLD) completes the 20-second cycle. At 20 seconds, a prolonged "beep" sound will indicate cycle completion. The Dual Cycle modality combined with IN-Motion™ technology offers more comfortable, effective and safer treatments.

## Accent<sup>XL</sup> Specifications

<b>Electrical</b>	Single phase, 208-240 Vac or 100-120 Vac Frequency of AC-line: 50/60 Hz
<b>System Dimensions</b>	(L × W × H): 55 × 40 × 100 cm
<b>Weight</b>	~ 50 kg



[www.almalasers.com](http://www.almalasers.com)

©2008 Alma Lasers, Ltd. All rights reserved. Alma Lasers, Ltd., its logo, Accent<sup>XL</sup>, IN-Motion™ and UniPolar™ are trademarks of Alma Lasers, Ltd. Product specifications are subject to change without notice.

### United States

485 Half Day Road, Suite 100  
 Buffalo Grove, IL 60089  
**Phn** (224) 377-2000  
**Fax** (224) 377-2050  
[contact@almalasers.com](mailto:contact@almalasers.com)

### Headquarters

14 Halamish Street  
 Caesarea Industrial Park  
 Caesarea, Israel 38900  
**Phn** +972-4-627-5357  
**Fax** +972-4-627-5368